



Disposal

The power source of this product is nickel-metal hydride battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. If you have recycling questions, contact your local waste management office.

Removing and recycling the battery

Do not disassemble the product at any time other than when disposing of it. Doing so earlier will destroy the product. Observe basic safety precautions. The appliance must be disconnected from the supply mains when removing the battery.

- 1 Before removing the battery, be sure to completely drain the unit of power.
- 2 Use a screwdriver to remove the screw cover and battery door.
- 3 Hold the handle upside down, press the toothbrush shaft, and remove the two hooks to release the internal components of the handle.
- 4 Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.
- 5 Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

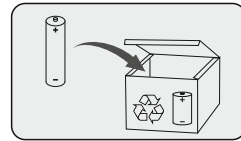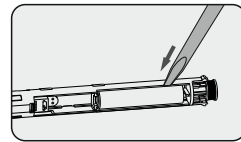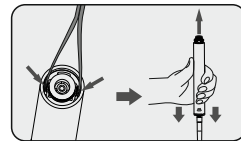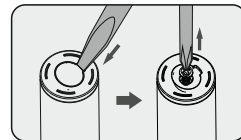Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.

Dr. Mayer
Life and Health

E-mail: office@doctor-mayer.com
Web: www.doctor-mayer.com


Dr. Mayer
LIFE & HEALTH

GTS2020

Sonic Toothbrush

Operating Instructions

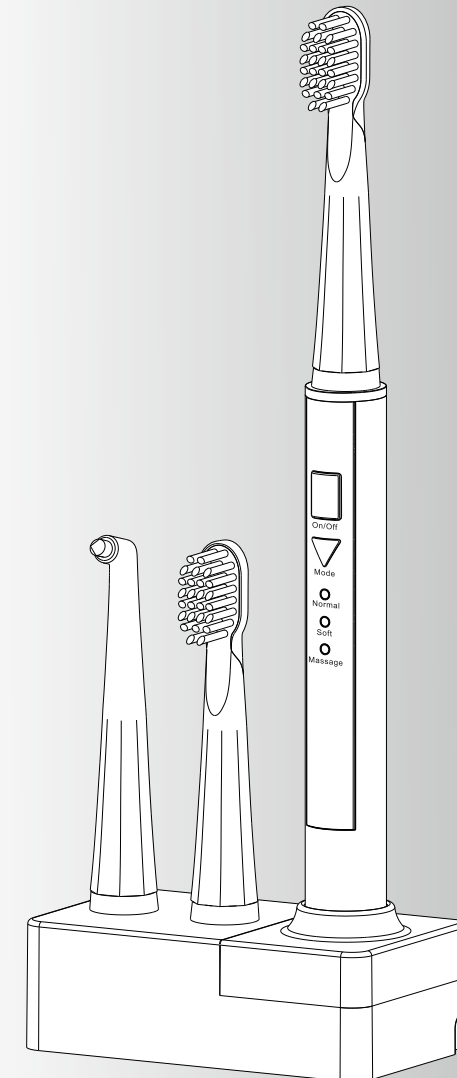ENG
RO
BG
HU

Before operating this product, please read these instructions completely and save this manual for future use.

Safety Instructions

Please read the safety instructions before using this product. The safety instructions provided here are to help you use this product safely and correctly, and also to prevent you from injuring yourself or others.

Dangers-Water

Keep the charger away from water:

- 1 Do not use the charger in the bath.
- 2 Do not immerse the charger into water or wash it with water or other liquid.
- 3 Do not place or store the charger where it can fall or be pulled into water.
- 4 Do not reach for the charger that has fallen into water. Unplug immediately to avoid the risk of electrocution.
- 5 Do not touch the plug with wet hands.

Dangers-Voltage

Check if the voltage indicated on the bottom of the charger corresponds to the local mains voltage before you connect the appliance:


- 1 The supply cord of this appliance cannot be replaced; if the cord is damaged, the appliance must be discarded. Do not use this product if it is not working properly, if it has been dropped or damaged, or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- 2 Do not damage, break or modify the cord. Do not bend, pull or twist the cord with unnecessary force. Do not place heavy objects on the cord or allow the cord to get pinched between objects.
- 3 Always have the charger replaced with one of the original type in order to avoid a hazard. Voltage converters do not guarantee compatibility.
- 4 Unplug the charger before performing any maintenance.
- 5 Do not use the charger outdoors or near heated surfaces.

Warnings

- 1 Use this product only for its intended use as described in this manual. Only use accessories that are recommended by the manufacturer.
- 2 Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- 3 Children shall not play with the product.
- 4 The toothbrush is only to be used with the charger provided with the product.
- 5 Do not attempt to remove the battery.
- 6 Do not place any metal objects, such as coins and clips, on the charger.
- 7 If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid plastic cracking.

- 8 Do not clean the toothbrush in the dishwasher.
- 9 The toothbrush is a personal care product and is not intended for use on multiple patients in a dental office or institution.
- 10 If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- 11 If you have had oral surgery or gum surgery within the past 2 months, consult your dental professional before using this product.
- 12 The toothbrush has been tested and is compliant with safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.

Configuration

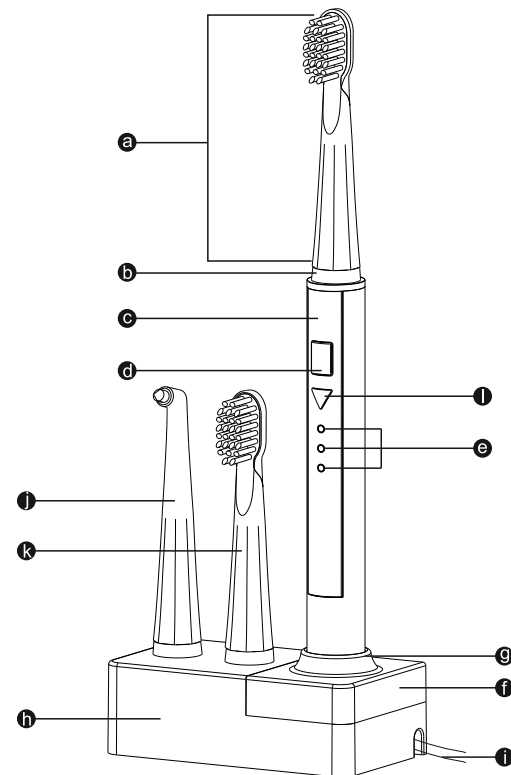
- | | |
|------------------------|---------------------|
| a. Brush head | g. Charging section |
| b. Identification ring | h. Charger |
| c. Toothbrush Handle | i. Power cord |
| d. On/Off button | j. Slot brush |
| e. Charging indicator | k. Super brush head |
| f. Charger cover | l. Mode switch |

Getting Started


Installing the charger Power Cord

Place the cord in the right, left or center cord guide, depending on where the charger will be placed.

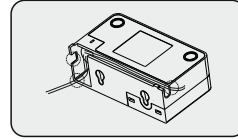
Installing near a sink or on a shelf

Place the charger on a level, secure surface. Insert the plug into an electrical outlet.

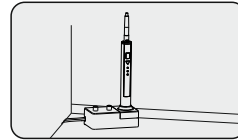
Mounting on the wall

Fasten two round wood screws on the wall, so that they extend about 2.5mm from the wall and are positioned 41mm apart.

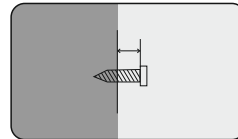
Align the holes on the back side of the charger with the screws and secure in place.

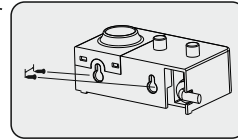Charging the toothbrush

Be sure to turn off the toothbrush before charging.

Place the handle in the charger. The charging indicator will light up green.

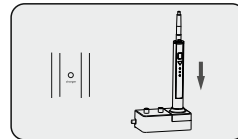
Notes:


- 1 It takes at least 24 hours to charge the battery fully. The indicator light does not go off even after fully charged.
- 2 When fully charged, the toothbrush can be used at least 40 minutes. The operating duration will become shorter as the battery gets older.
- 3 The toothbrush must be charged for 24 hours, if:
 - a. this is the first time using the toothbrush;
 - b. the toothbrush has not been used for more than three months; or
 - c. if the charging indicator blinks green during use.
- 4 To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.
- 5 During charging, the toothbrush will get warm, but this is normal.
- 6 If the toothbrush operates only for 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle. Take the unit to an authorised service center to have the battery replaced.

How To Use

- 1 Align the slot of brush head to the protrusion of handle and slide down, and then turn the brush head until it clicks into place.

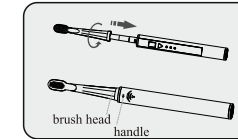- 2 Wet the bristles and apply a small amount of toothpaste.

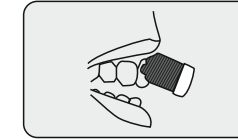
- 3 Place the bristles against the teeth at a slight angle towards the gum line.

- 4 Press the mode selector to the mode.
 - a. Normal: Standard mode for superior teeth cleaning.
 - b. Soft: Gently cleans your teeth or any other sensitive areas.
 - c. Massage: Pulsates for gentle gum stimulation.


- 5 Press the on/off button to turn on the toothbrush.

- 6 Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.


- 7 In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use.

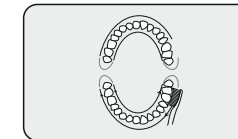- 8 Turn the brush head in direction as shown, and then pull it to detach from the toothbrush.

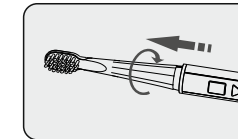- 9 Rinse both the brush head and handle under running water wipe them dry.

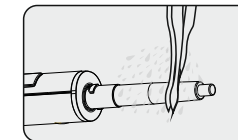
- 10 Put the handle back into the charger.

Maintenance

Be sure to remove the plug from the electrical outlet before cleaning the product.

Charger

- 1 Pull the charger cover to remove.
- 2 Using a soft cloth, wipe off any dirt on the charger and the charger cover.
- 3 Clean the plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheat or fire.

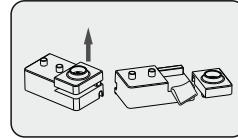Toothbrush

- 1 Remove the brush head.
- 2 Wipe all surface areas with a damp cloth and place the handle back into the charger. Do not submerge the handle into water.
- 3 Rinse the brush head and bristles after each use.
- 4 Do not clean in the dishwasher.

Troubleshooting

Problem	Possible cause	Solution
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 24 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Take the unit to an authorised service center to have the battery replaced.
	Charging time was insufficient.	Charge the main unit for at least 24 hours.

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